



SYSTEM OF GP CONTACT LENS WITH  
ECCENTRICITY CONTROL AND  
PERIPHERAL SUPPORT

S P E C T R U M  
I N T E R N A T I O N A L

Myopia - Hyperopia  
KC - Irregular Corneas  
Astigmatism - Intolerance  
DMP - Presbyopia



Made by:



# 1. Introduction

## Basic Concepts



- Corneal Physiology and form coefficients
- Design of the Elipsys lens
- Material
- Data Sheet

## 2. How to Fit Elipsys

### Types of fitting

- **Regular cornea**
  - Initial lens calculation
  - Assessment of fitting
- **Irregular Cornea**
  - Initial lens calculation
  - Assessment of fitting

# 3. Elypsis Multifocal

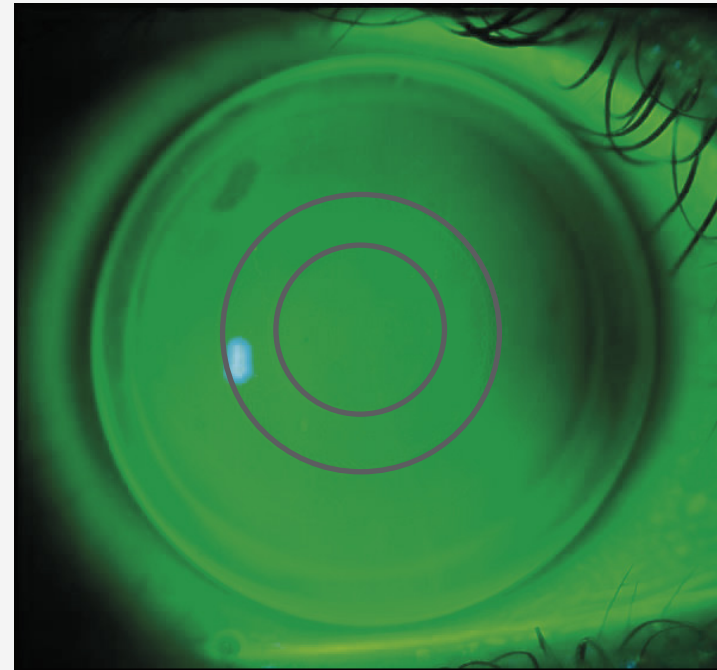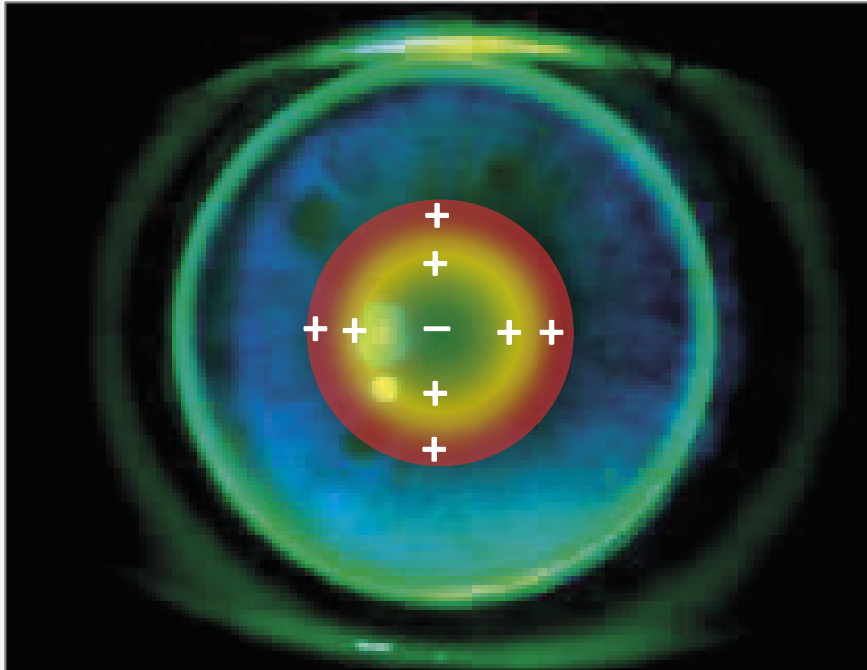


# 4. Introduction: Basic Concepts


## Corneal Physiology and form coefficients

- P Value (p)
- Asphericity (Q)
- Eccentricity (e)
- Shape Factor (SF)

	p	Q	e	SF
p	-	$1+Q$	$1-e^2$	$1-SF$
Q	$p-1$	-	$-e^2$	$-SF$
e	$\sqrt{1-p}$	$\sqrt{-Q}$	-	$\sqrt{SF}$
SF	$1-p$	$-Q$	$e^2$	-

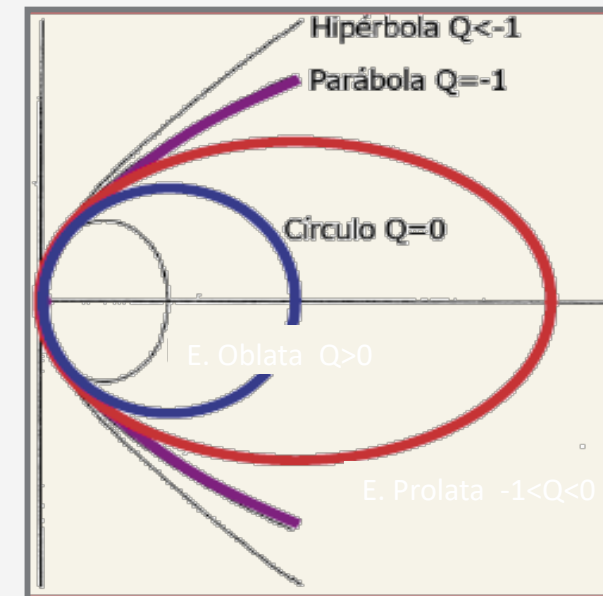

# Design of the Elipsys lens

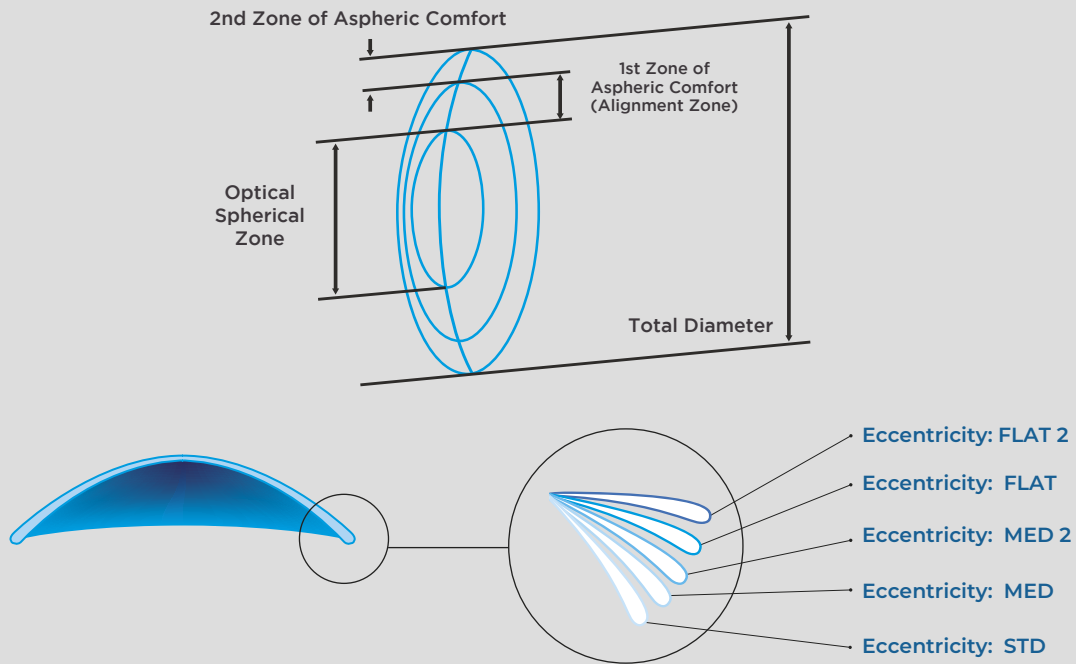

Sphero-Aspheric design on the posterior section with peripheral support control that allows a better corneal alignment optimizing the tear exchange, maintaining its centering and stability.

Multisphere geometry on the anterior side.

Indicated for patients with regular corneas or irregular corneas with keratoconus or other corneal ectasias.

Design with optical zone of variable diameter that allows it to stay centered and stable in most irregular corneas.

The new concept of “peripheral support control” is ideal for corneas that have different levels of eccentricity. Based on a statistical analysis, a nomogram was created identifying several levels of eccentricity for the different corneal conditions.

# Material

**Contamac**<sup>™</sup>  
EMPOWERING YOUR VISION

**OPTIMUM**  
EXTRA

Property	Optimum Extra
Classification	Focon III 4
USAN	Roflufocon D
Oxygen Permeability (Dk)	100 barrers
Refractive index	1,431
Hardness (Shore D)	75
Retracting contact angle	3°
Visibility Tint	Blue
Recommended Replacement	18 - 24 months

# Data sheet

## DATA SHEET

Geometry (interior):	Spherical-Aspheric with peripheral alignment control
Geometry (exterior):	Sphero-progressive

### Characteristics:

### Available Parameters:

Power (D):	25D
Eccentricity (Ecc):	STD, MED, MED 2, FLAT and FLAT 2
Base curve (mm):	5.00mm to 11mm (0.05mm steps)
Total diameter (mm):	7.00mm to 12.60mm (0.10mm steps)
Optic zone diameter (mm):	6.00mm to 10.00mm (0.10mm steps)
Central thickness:	0.17 (-3.00D) standard
Material:	Optimum Extra - CONTACMAC
Classification:	Focon III 4
Permeability to oxygen:	DK 100
Visibility Tint:	Blue
Recommended replacement:	Between 18 and 24 months
Plasma Treatment:	Yes
Identification Engraving:	Yes
Presentation:	Vial with security seal
Maintenance System:	Preservative + Isopropyl alcohol cleaner

### Fitting Set with 30 Lenses

ElipSYS STD	Regular Cornea	DT. 9.80mm 10 Lenses
ElipSYS KC	Irregular Cornea	DT. 9.60mm 10 Lenses
ElipSYS SE	Readjustments / Decentralizations	DT. 11.20mm 10 Lenses



## 2. How to Fit Elipsys

### Important:


#### Before starting, make sure you have:

- Update graduation, AVsc and AVcc
- Current Keratometry or topography
- Horizontal diameter of visible iris (DHIV)
- Pupillary diameter (adap. subpalpebrales, LC MF, etc..)
- Evaluation of the patient's tear
- Valuation of the corneal surface (Avoid surprises)
- A good anamnesis (previous adaptations, expectations, etc.)

**KNOW YOUR PATIENT**


# Subpalpebral Fitting

- Ø Large posterior optical zone
- Wide bands
- Base radio equal to and / or greater than Kflat

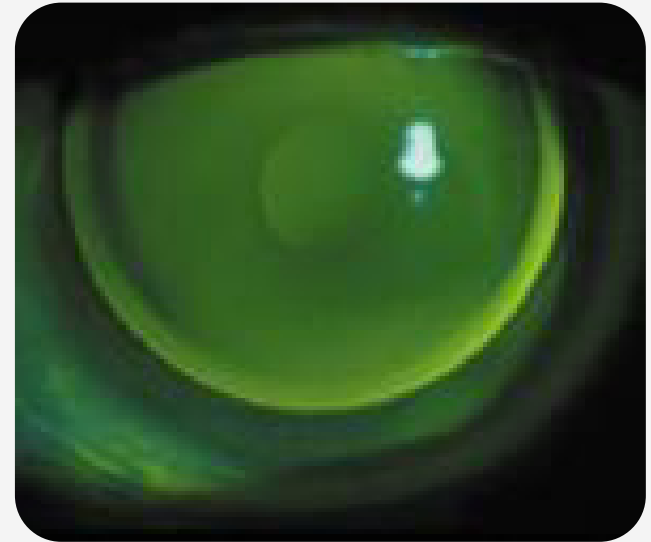

# Subpalpebral Fitting

- $\emptyset$  Large posterior optical zone
- Narrow peripheral bands
- Base radio equal to or less than Kflat
- $\emptyset$  total < Parpebral opening

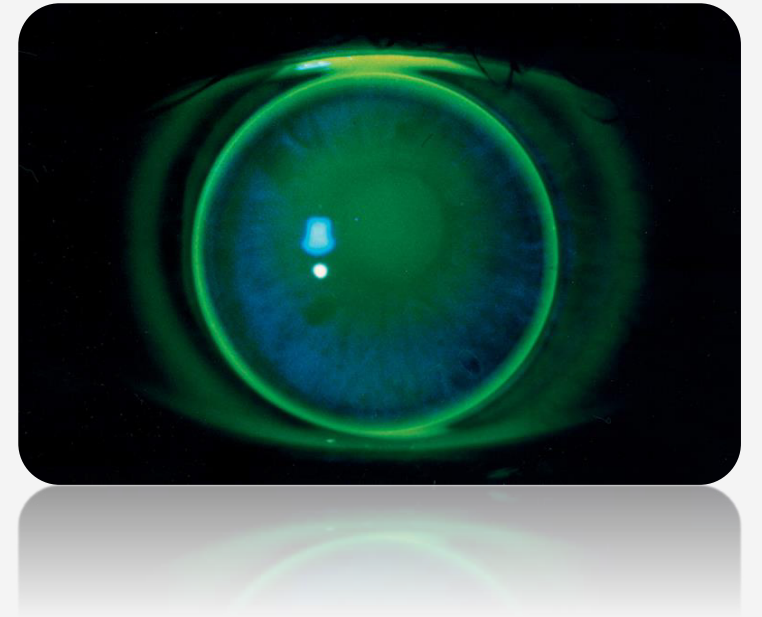

# Córnea Regular and Irregular

## *Elypsis Family (STD, KC and SE)*


**Step N° 1:** *Calculate the total diameter of the lens*



9.60 mm  
Diameter



9.80 mm  
Diameter

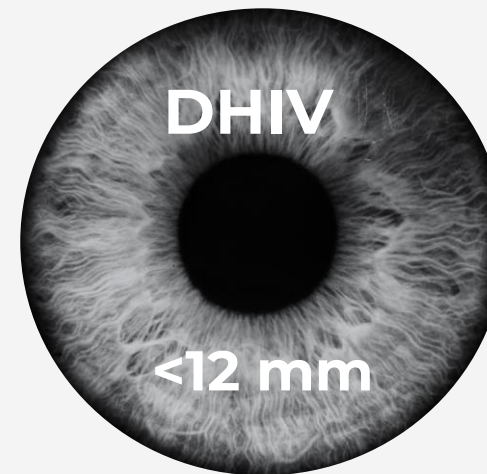

Between 9.80 mm and  
DHIV - 14.0mm

# Regular Cornea

## *Elipsys STD and SE*

### Step N° 2: Selection of the base curve

Select the lens from Fitting Set, with the base curve that is closer to the value obtained according to the calculations indicated by the table:

Elipsys STD and Elipsys SE			
AC: 0.00 to 0.75 D		AC: 0.75 to 1.75 D	
Diameter LC	BC	Diameter LC	BC
9.60	$K_{\text{flat}} + 0.05 \text{ mm}$	9.60	$K_{\text{flat}}$
9.80	$K_{\text{flat}} + 0.05 \text{ mm}$	9.80	$K_{\text{flat}}$
11.20	$K_{\text{flat}} + 0.25 \text{ mm}$	11.20	$K_{\text{flat}} + 0.20 \text{ mm}$

Design NOT recommended for values  $> 2.00 \text{ D}$  dof corneal astigmatism (regular) It would condition the focus, comfort and visual stability.

# Regular Cornea

## *Elypsis STD and SE*

### Step N° 2: Selection of Base Curve

- For values > 2.00 D of corneal astigmatism (regular) or astigmatism 2.00 D but with an extension of limbus to limbus.

### Treat the case as if it was an Irregular Cornea

- Select the lens out of the fitting set, with the base curve that is closer to the value obtained according to the calculations indicated by the table:

Diameter	BC
9.60 y 9.80	$\frac{K_{\text{flat}} + K_{\text{steep}}}{2} - 0,10$
11.20	$\frac{K_{\text{flat}} + K_{\text{steep}}}{2} + 0,10$

### Analyze fluorogram and evaluate:

Centering and lifting of the edge

- Low edge lifting → ask for Exc. “MED”
- Excessive edge lifting → Try a steeper BC

# Regular Cornea

## *Elypsis STD and SE*

### Step N° 1



Between 9.80 mm and  
DHIV - 14.0mm

### Step N° 2

42.19 (8.00) @ 180 / 42.92 (7.86)@ 90  
 $K_{\text{flat}} = 42.19 \text{ D (8.00 mm)}$   
Corneal Cyl= 0.73 D  
Base Curve=  $K_{\text{flat}} + 0.05 \text{ mm}$   
= 8.05 mm  
Base Curve=  $K_{\text{flat}} + 0.25 \text{ mm}$   
= 8.25 mm

**First Fitting Lens:** 8.00/-0.50/9.80/STD (*Design Elypsys STD*)  
8.30/+1.00/11.20/STD (*Design Elypsys SE*)

# Regular Cornea

## *Elypsis STD and SE*

### **Step N° 3:** *Evaluate centering and fluorogram*

**1. Static Behavior** (holding eyelids and focusing the lens):

- Movement of the contact lens without eyelid action (Push Up)
- Static fluorogram (acceptable, open, closed)

**2. Dynamic Behavior** (under the action of the eyelids):

- Movement
- Centering
  - Increase diameter to raise the lens
  - Decrease diameter to lower the lens

If the GP contact lens provides good AV, comfort and absence of problems, the fitting is good




# Regular Cornea

## *Elypsis STD and SE*



**Ideal**

- ✓ Uniform distribution of fluorescein
- ✓ No touch zones
- ✓ Right edge lifting
- ✓ Good movement and centered

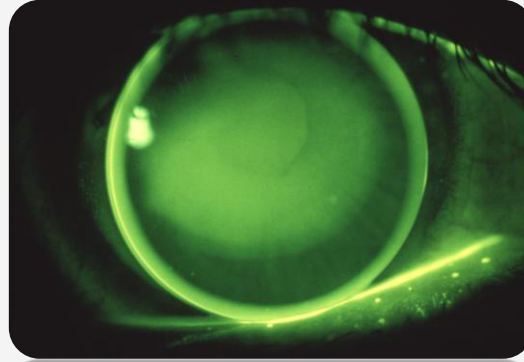

**Steep**


### STEP LENS

Apply the solutions in this order

Make BC Flatter

Decrease Diameter

Decrease Optic DZ

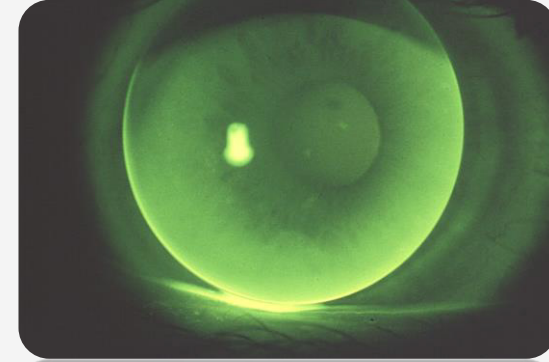

**Flat**

### FLAT LENS

Apply the solutions in this order

Make BC Steeper

Increase Diameter

Increase Optic DZ

# Irregular Cornea

## *Elypsis KC and SE*

### Step N° 2: Selection of Base Curve

Select lens out of the fitting set, with the base curve that is closer to the value obtained according to the calculations indicated by the table:

ElipSYS KC and ElipSYS SE	
Diameter	BC
9.60 and 9.80	$\frac{K_{\text{flat}} + K_{\text{steep}}}{2} - 0,10$
11.20	$\frac{K_{\text{flat}} + K_{\text{steep}}}{2} + 0,10$

# Irregular Cornea

## *Elypsis KC and SE*

### Step N° 1



Between 9.80 mm and  
DHIV - 14.0mm

### Step N° 2

49.63 (6.80) @ 180 / 53.57 (6.30)@ 90

$$\frac{K_{\text{flat}} + K_{\text{steep}}}{2} - 0,10 = 6.45\text{mm}$$

$$\frac{K_{\text{flat}} + K_{\text{steep}}}{2} + 0,10 = 6.65\text{mm}$$

**1st fitting lens:** 6.40/-10.75/9.60/FLAT (*Design Elypsys KC*)  
6.70/-8.75/11.20/FLAT 2 (*Design Elypsys SE*)

# Irregular Cornea

## *Elypsis KC and SE*

### **Step N° 3:** *Evaluate centering and fluorogram*

#### **1. Static Behavior** (holding eyelids and focusing the lens):

- Movement of the contact lens without eyelid action (Push Up)
- Static fluorogram (acceptable, open, closed)

#### **2. Dynamic Behavior** (under the action of the eyelids):

- Movement
- Centering
  - Increase diameter to raise the lens
  - Decrease diameter to lower the lens

If the GP contact lens provides good AV, comfort and absence of problems, the fitting is good

# Irregular Cornea

## *Elypsis KC and SE*

### **Step N° 3:** *Evaluate centering and fluorogram*


- **IMPORTANT**

- Initially, evaluate the center of the lens
  - Make the BC steeper or flatter until getting the ideal fitting
- Evaluate the lifting of the edge
  - Adjust with the eccentricity (see troubleshooting).

# Irregular Cornea


## *Elypsis KC and SE*

### Ideal

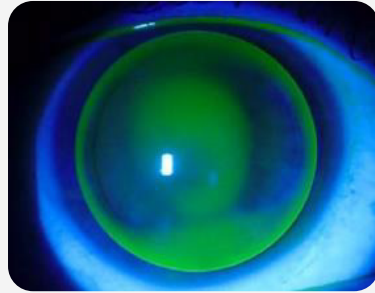

- ✓ Apical light touch (toque pluma)
- ✓ Middle-periphery support
- ✓ Right edge lifting
- ✓ Good movement and centered

### Steep

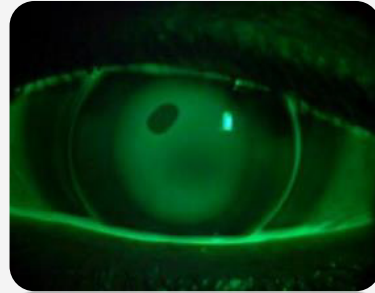

#### STEP LENS

Apply the solutions in the order

Make BC Flatter

Increase Eccentricity

Decrease Diameter

Decrease Optic DZ

### Flat



#### FLAT LENS

Apply the solutions in the order

Make BC Steeper

Decrease Eccentricity


Increase Diameter

Increase Optic DZ

# Irregular Cornea

## *Elypsis KC and SE*

*Ideal adaptation technique in irregular cornea*

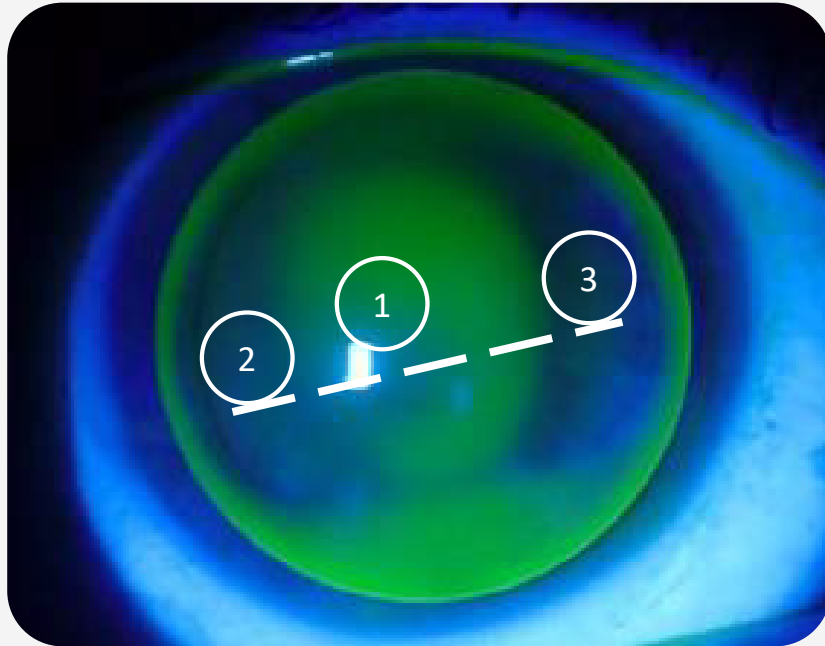

### Three Support Points

- 1 Slight Apical Touch (feather touch)
- 2 y 3 Support in the Mid-periphery


# Irregular Cornea

## *Elypsis KC and SE*

### Troubleshooting

If you observe:

Lifting of the Edge	Solution
<b>Insufficient</b> (less than de mm)	<b>Increase Eccentricity</b> (ei.: MED → MED 2, FLAT or FLAT 2)
<b>Excessive</b> (greater than 1 mm)	<b>Decrease Eccentricity</b> (ei.: FLAT 2 → FLAT, MED 2 or MED)

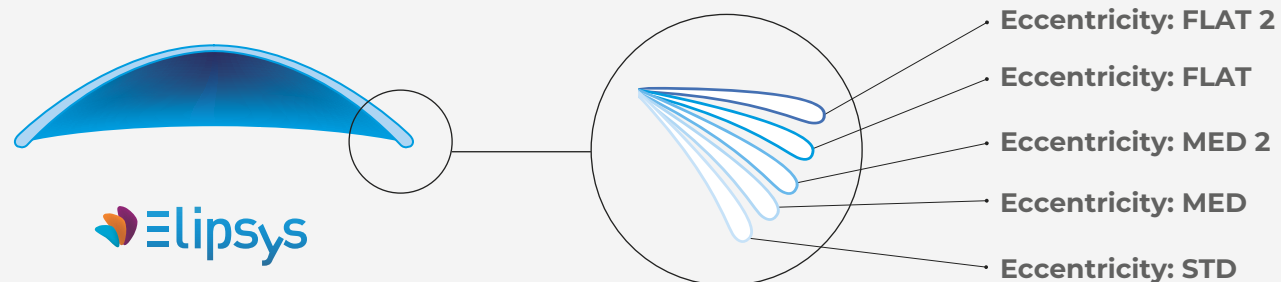



# Regular and Irregular Cornea

## *Elypsis Family STD, KC and SE*

### Troubleshooting

	Base Curve	Diameter	Eccentricity
Dropped Lens	Flatten	Increase	Increase
Lens riding High	Close	Decrease	Decrease
Central Staining	Close		
Staining 3 and 9		Decrease	Increase
Superior Limb Stain		Decrease	
Blurry Vision	Value centered and / or Pupular Diameter		
Poor Vision	Flatten		
Central Pooling	Flatten	Decrease	
Dimple Veiling (trapped bubbles)	Flatten	Decrease	

# Regular and Irregular Cornea

## *Elypsis Family STD, KC and SE*



### Step N° 4: *Power of the lens*

- Perform an Over-refraction
  - Do not forget to apply distometry (vertex distance) for values higher than  $\pm 3.75$ .
- Keep in mind that this design DOES NOT compensate for internal astigmatism
  - Ast. Rx= Corneal Astigmatism + Internal Astigmatism
- Apply rules:
  - Make Steeper add negative
  - Make flatter add positive


# Regular and Irregular Cornea

## *Elypsis Family STD, KC and SE*

### Step N° 4: *Power of the Lens*

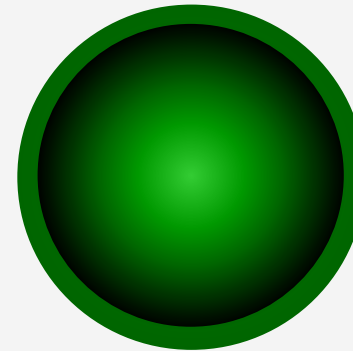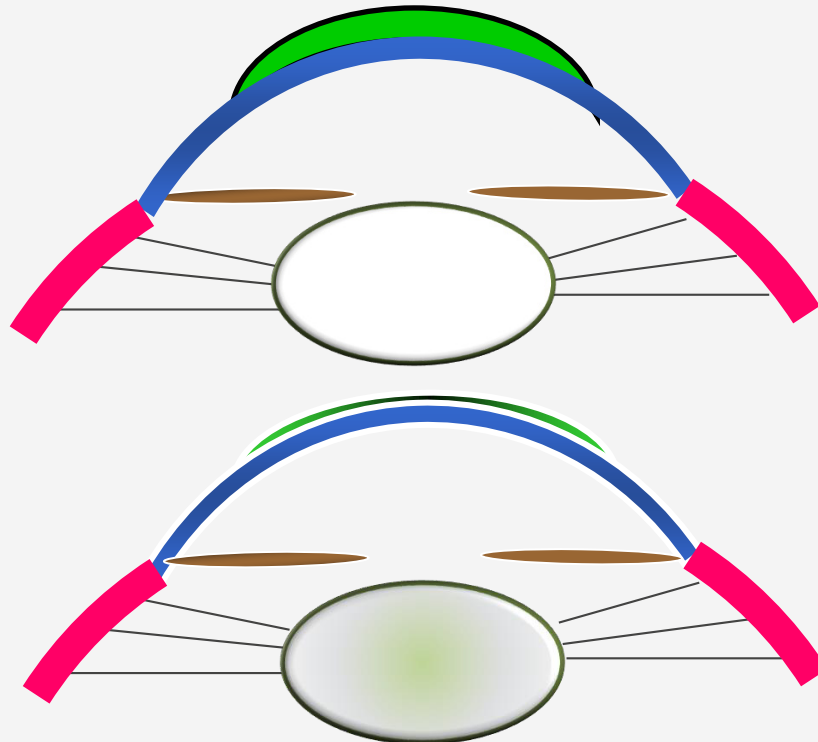

Make BC steeper  
Meniscus lacrimal positive  
Add Negative

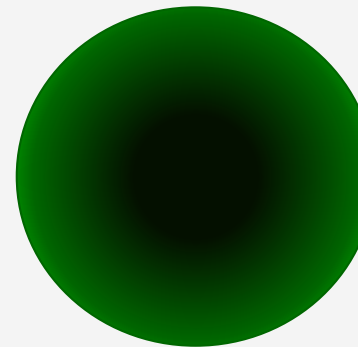

Make BC flatter  
Meniscus lacrimal negative  
Add Positive

# Regular and Irregular Cornea

## *Elypsis Family STD, KC and SE*

### **Step N° 5:** *Order of final lens (parameters to be identified)*

- Design/Geometry
  - Elipsys STD
  - Elipsys KC
  - Elipsys SE
- Base curve
- Power
  - Fitting lens number + Over refraction (Lab.will apply distometry)
  - Definitive lens power (distometry already applied)
- Diameter
- Eccentricity


# Regular and Irregular Cornea

## *Elypsis Family STD, KC and SE*

### Possibility of Multifocal Design

- Presbyopic patients wearing Elypsys family contact lenses can now enjoy clear vision at all distances, regardless of their corneal condition\*\* (regular and irregular cornea).
- Simultaneous vision multifocal design + Translation
  - Center - Far (CD)
- Required parameters for you order (empirically):
  - In addition the parameters required to order an Elypsys lens:
    - Pupillary Diameter
    - ADD
    - Ocular Dominance
    - (Sensorial by fogging test with + 2.00 D lens)

(\*\*) We recommend special attention when trying to pass, the cases of very irregular cornea, to a Multifocal geometry, since an accumulation of optical aberrations both corneal and the same multifocality of the lens, may condition the final visual result.

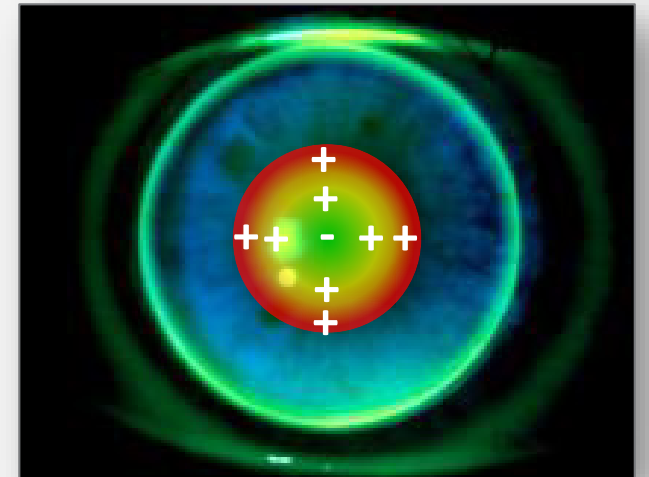


# Regular and Irregular Cornea

## *Elypsis Family STD, KC and SE*

### Multifocal design possible

- Possibility of fitting set with multifocal geometry.

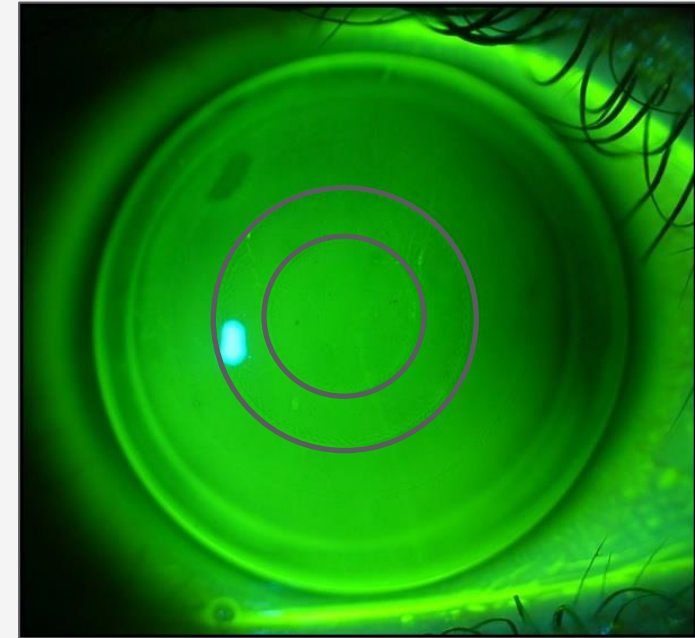



S P E C T R U M  
I N T E R N A T I O N A L

[www.spctinternational.com](http://www.spctinternational.com)

